

YOUR EXPERT IN DENTAL INSTRUMENTS



OTTO **LEIBINGER**
L-DENT

THE COMPANY

Otto Leibinger GmbH was established in 1978. Its founder and name-sake comes from a family with a long tradition in instrument manufacture. Our location in Mühlheim near Tuttlingen provides the best conditions for manufacturing high-quality dental hand instruments as the world's centre for medical technology. In collaborating with us, customers greatly appreciate our flexibility and service strengths. With direct, personal care and a solution-oriented way of thinking and treating, we set the course for a long-term, cooperative partnership.



It is our pleasure to be at your service for all your questions and special requests.

A handwritten signature in blue ink, appearing to read 'O. Leibinger'.

Otto Leibinger
President

A handwritten signature in blue ink, appearing to read 'Tina Leibinger-Toth'.

Tina Leibinger-Toth
Vice-President

EXPLORERS, FILLING INSTRUMENTS

When the company started 30 years ago, this product range has been the first range being manufactured as OEM for well known suppliers. Over the years our program has enlarged impressively and we are proud to say that we are one of the most experienced manufacturer in this field, supplying superior quality.



EXPLORERS & FILLING INSTRUMENTS

GRACEY CURETTES AND SCALERS


The whole range of Gracey Curettes for periodontal treatment, available with standard handle or handle with area-corresponding color coding, that helps you to quickly identify the area of application. Tips are being made of a special stainless steel that ensures a long lifetime and a long lasting sharpness.



SCALERS

2

GRACEY CURETTES

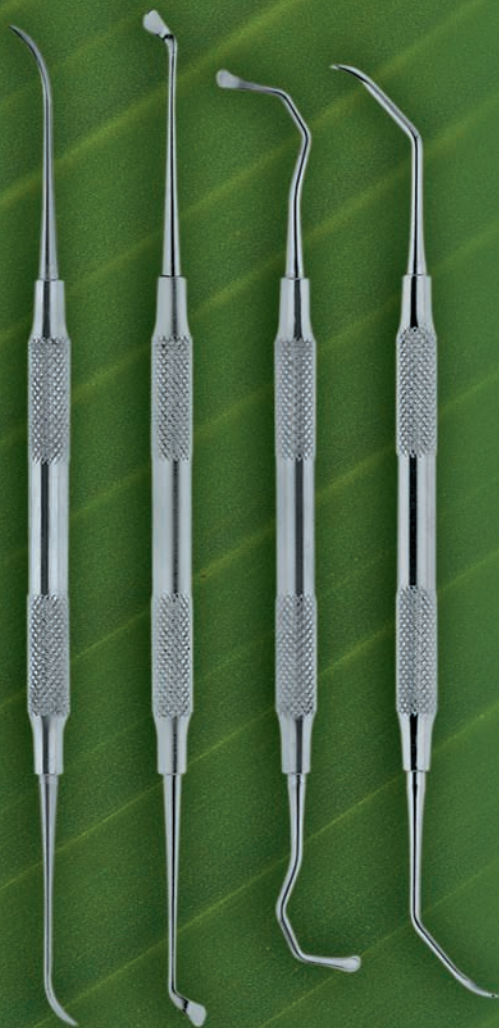

SINUS LIFT CURETTES

Available as set of 15 different pattern in rigid version and as a set of 6 different pattern in a flexible version.

OSTEOTOMES / CHISELS

In case of poor bone thickness up to 3 mm, it is possible to split the jaw-bone (Bone Split) and separate it by means of special chisels, the so called Osteotomes. Finally a gap occurs between buccal and lingual respectively buccal and palatal bone-lamella, where you can insert an implant. There are two sizes in straight and curved available.



CHISELS




OSTEOTOMES

OSTEOTOMES WITH ERGONOMIC HANDLE

For lateral bone condensation of the closed sinus elevation technique. Available in straight and bayonet shaped, with or without stop and adjusting screw. Tips come in pointed, convex or concave version.



OSTEOTOMES

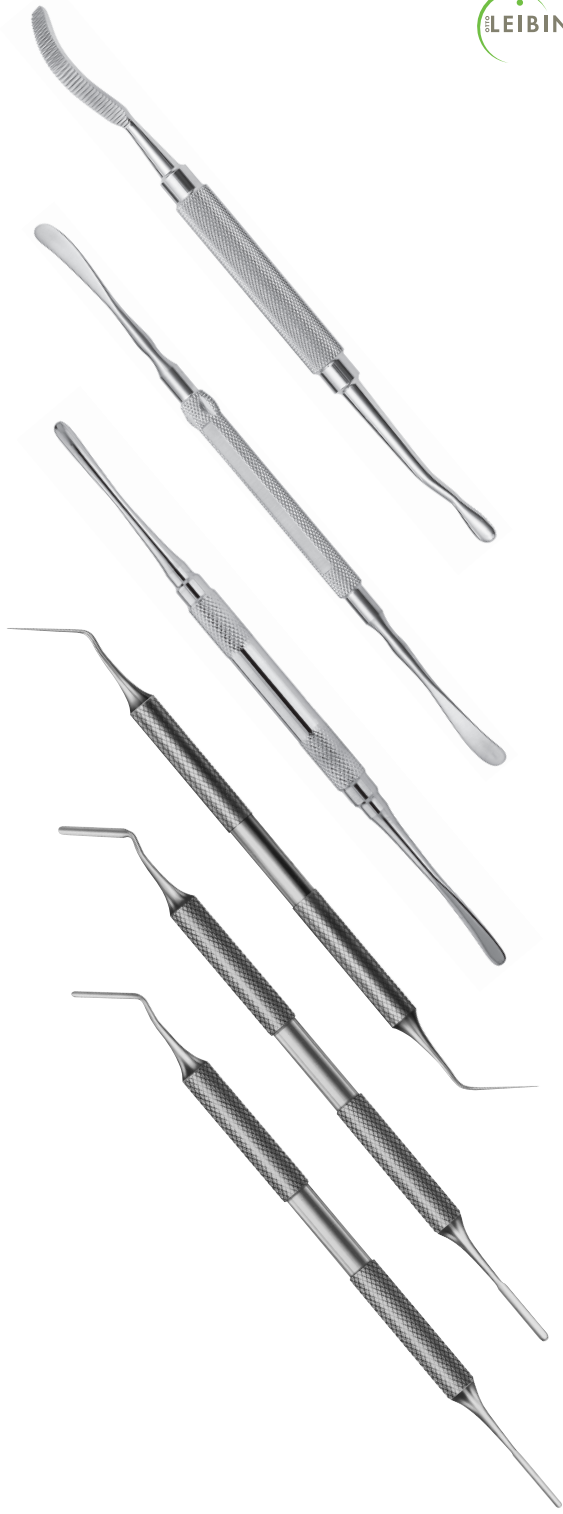

ELEVATORS & RASPATORIES

4

PERIOTOMES & BONE FILES

BONE MILL

Manual Bone Mill to chop bone fragments being removed from the oral cavity for use in regenerative bone surgery. This bone mill is consisting of a base cup that is made of stainless steel to ensure safe positioning on the working place, additionally ensured by a special anti-slip rubber base. The pestle and the bone collecting cup are made of titanium.



MEASURING GAUGE

Due to the multifunctional character it is easy to set an implant. The measuring rods allow a fast and defined control of parallel and symmetrical distances between implants in distal and mesial direction. The hole at the end of the sliding callipers serves as a drilling gauge, drilling guide and depth measurement.





STEAMER NET RICHTER

5

STEAMER NET RICHTER

An innovative product to be used in dental labs. Essential for preventive and safe steam cleaning of small workpieces. Developed in cooperation with an experienced user.

DENTAL HF-ELECTRODES
FOR CUTTING
AND COAGULATING

Suitable for all handles with 1,6mm shank holder,
available in rigid and flexible version. (different
color of isolation)



HF-ELECTRODES

THE RANGE OF PRODUCTS

Beside all these specialities that are being manufactured in our facility in Muehlheim/Germany, we also offer a complete program of all kinds of dental hand instruments, of course all of excellent quality and competitive pricing as there are:

- EXTRACTING FORCEPS AND ROOT ELEVATORS
- ORTHODONTIC PLIERS
- FULL RANGE OF INSTRUMENTS AND TOOLS FOR THE DENTAL LAB/DENTAL TECHNICIAN
- SCISSORS, NEEDLE HOLDERS, TWEEZERS, RETRACTORS AND A WIDE VARIETY OF INSTRUMENTS FOR DENTAL SURGERY AND MANY MORE
- OEM PRODUCTION OF STANDARD AND CUSTOMIZED ITEMS



OTTO LEIBINGER GMBH
GRIESWEG 27
78570 MÜHLHEIM / GERMANY

FON: +49 (0) 7463-7232
FAX: +49 (0) 7463-222

info@otto-leibinger.de
www.otto-leibinger.de